**EN**

## **DeVilbiss® Pulmo-Aide® Compact Compressor/Nebulizer Instruction Guide**

**Models 3655D and 3655I**

**CAUTION**—Federal (U.S.A.) law restricts this device to sale by, or on the order of a physician.

Assembled in the USA. Not made with natural rubber latex.

Models 3655D Intertek ETL Certified to CAN/CSA C22.2 No. 601.1 and UL60601-1.

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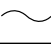
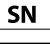
**IEC SYMBOLS**

	Attention - Consult Instruction Guide		Type BF equipment– applied part
	Consult instructions for use		Double Insulated
	Alternating Current		Serial Number
	This device contains electrical and/or electronic equipment that must be recycled per EU Directive 2012/19/ EU-Waste Electrical & Electronic Equipment		

**IMPORTANT SAFEGUARDS**

When using electrical products, especially when children are present, basic safety precautions should always be followed. Read all instructions before using. Important information is highlighted by these terms:

**DANGER**– Urgent safety information for hazards that will cause serious injury or death.

**WARNING**–Important safety information for hazards that might cause serious injury.

**CAUTION**– Information for preventing damage to the product.

**NOTE**– Information to which you should pay special attention.

**READ ALL INSTRUCTIONS BEFORE USING.**

**DANGER**

To reduce the risk of electrocution:

1. Always unplug this product immediately after using.
2. Do not use while bathing.
3. Do not place or store product where it can fall or be pulled into a tub or sink.
4. Do not place in or drop into water or other liquid.
5. Do not reach for a product that has fallen into water. Unplug immediately.

## **WARNING**

**To reduce the risk of burns, electrocution, fire or injury to persons:**

1. The product should never be left unattended when plugged in.
2. Close supervision is necessary when this product is used by, on, or near children or physically incapacitated individuals.
3. Use this product only for its intended use as described in this guide. Do not use attachments not recommended by the manufacturer.
4. Never operate this product if:
  - a. It has a damaged power cord or plug.
  - b. It is not working properly.
  - c. It has been dropped or damaged.
  - d. It has been dropped into water.Return the product to an authorized DeVilbiss service center for examination and repair.
5. Keep the power cord away from heated surfaces.
6. Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
7. Never use while drowsy or asleep.
8. Never drop or insert any object into any opening or hose.
9. Do not use outdoors. This product is for household use only.
10. Do not use in an oxygen enriched environment.
11. Do not operate unit in a backpack or carrying bag.
12. Care should be taken when this product is used in the vicinity of other electronic equipment concerning the possible interference between them. They may need to be moved apart.

**NOTE**– *This unit employs a temperature and/or current limiting device that will shut off the entire product should a malfunction occur or in the event of product misuse, abuse, or abnormal operation. Contact your qualified DeVilbiss provider for a replacement unit if necessary.*

This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature. If repair or replacement of the cord or plug is necessary, contact your qualified DeVilbiss provider.

**NOTE**– *The exterior surfaces of the unit and the air outlet connector may become slightly warm to the touch during use. Warm air can also be felt exiting the vents on the sides of the case of the unit. This is normal.*

## **INTRODUCTION**

Your doctor has prescribed a liquid medication to treat your respiratory condition. To best utilize this liquid medication, he has prescribed a DeVilbiss brand compressor/nebulizer. Your DeVilbiss compressor/nebulizer works to convert the medication to a high-quality mist of fine particles that penetrates deep into the lungs.

Be sure that you read and understand the information in this instruction guide. By following these simple instructions and the advice of your physician, your compressor will become an effective addition to your therapeutic routine.

### **Statement of Intended Use**

The DeVilbiss compressor/nebulizer includes an AC powered air compressor that provides a source of compressed air for home health care use. The compressor is used in conjunction with a jet (pneumatic) nebulizer to convert certain inhalable drugs into an aerosol form for inhalation by a patient. The device can be used with adult or pediatric patients.

**NOTE**– *If using a nebulizer other than one distributed by DeVilbiss, verify that the pressure and/or flow characteristics of the nebulizer (while driven by the DeVilbiss compressor) meet the nebulizer manufacturer's recommendation.*

## IMPORTANT PARTS

### DeVilbiss 3655 Compressor (Figures A & B)

1. Power switch
2. Filter (non-washable)
3. Air-outlet connector
4. Carrying handle
5. Power cord storage recess

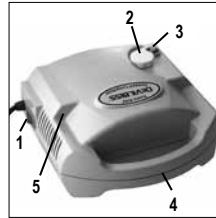
Figure A

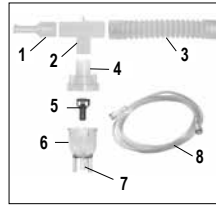
Figure B


### Disposable Nebulizer

1. Mouthpiece
2. T-piece
3. Corrugated Tubing
4. Cap
5. Baffle
6. Medication Cup
7. Nebulizer Air-Inlet Connector
8. Tubing

Nebulizer parts are not sold separately. Additional nebulizers (part #3655D-621) are available through your DeVilbiss provider.

Disposable Nebulizer

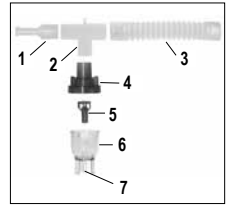Reusable Nebulizer

### Reusable Nebulizer (Optional)

1. Mouthpiece
2. T-piece
3. Corrugated Tubing
4. Cap
5. Baffle
6. Medication Cup
7. Nebulizer Air-Inlet Connector

Nebulizer parts are not sold separately. Additional nebulizers (part #900D-610) are available through your DeVilbiss provider.

## HOW TO OPERATE YOUR COMPRESSOR

**NOTE**— Before initial operation, your nebulizer should be cleaned following cleaning instructions, or as recommended by your physician or DeVilbiss provider.

1. Place compressor on a level, sturdy surface so that the controls can be easily reached when seated. Never block the air openings of the product or place it on a soft surface, such as a patient's lap or a bed or couch where the air openings may be blocked.
2. Make sure the power switch is in the "Off" position (Fig. 1). Unwrap power cord and plug power cord into an appropriate wall outlet (refer to Specifications). Plug model 3655 into a grounded wall outlet.

Figure 1

### DANGER

**DeVilbiss compressor/nebulizer must be operated on the specified power source to avoid risk of electric shock and damage to the compressor.**

3. Wash hands.
4. Connect one end of the nebulizer tubing to the compressor air-outlet connector (Fig. 2).

**NOTE**— During high-humidity weather conditions, condensation (water build up) can occur in the nebulizer tubing.

5. Assemble clean nebulizer and add prescribed medication (Fig. 3).
6. Assemble mouthpiece, t-piece, and corrugated tubing (if applicable) and insert into the top of the nebulizer cap. If using an aerosol mask, insert the bottom part of the mask directly into the top of the nebulizer cap.
7. Attach tubing to nebulizer air-inlet connector (Fig. 4).
8. Press the power switch "On" (Fig. 1) to start the compressor.

Figure 2

9. Begin treatment by placing the mouthpiece between teeth. With mouth closed, inhale deeply and slowly through mouth as aerosol begins to flow, then exhale slowly through the mouthpiece (Fig. 5). If treatment needs to be interrupted, simply press power switch "Off."

**NOTE**—Your physician may recommend a specific breathing pattern during treatment. If so, follow his/her recommendations.

10. If an aerosol mask is used, place the mask over mouth and nose (Fig. 6). As aerosol begins to flow, inhale deeply and slowly through mouth, then exhale slowly.
11. When treatment is complete, turn unit off by pressing the power switch to the "Off" position. Unplug the unit from the wall power outlet.

## NEBULIZER CLEANING

All parts of the nebulizer, except tubing, should be cleaned according to the following instructions. Your physician and/or DeVilbiss provider may specify a certain cleaning procedure. If so, follow their recommendations.

### WARNING

To prevent possible risk of infection from contaminated medication, cleaning of the nebulizer is recommended after each aerosol treatment. Disinfecting is recommended once a day.

### Clean After Every Use:

1. With power switch in the "Off" position, unplug power cord from wall outlet.
2. Disconnect tubing from the air-inlet connector and set aside.
3. Disassemble mouthpiece or mask from cap. Open nebulizer by turning cap counterclockwise and removing baffle.
4. Wash all items, except tubing, in a warm water/dishwashing detergent solution. Rinse under warm tap water for 30 seconds to remove detergent residue. Allow to air dry.
5. For reusable nebulizer only, clean daily in dishwasher using top shelf.

**CAUTION**—The reusable nebulizer is dishwasher safe, but do not place any disposable nebulizer parts in an automatic dishwasher; doing so could cause damage.

### Disinfect Daily:

1. Using a clean container or bowl, soak items in three parts warm water (131°F-149°F/55°C-65°C) to one part white vinegar (>=5% acetic acid concentration) for 30 minutes (Fig. 7) OR use a medical bacterial-germicidal disinfectant available through your provider. Be sure to follow manufacturer's instructions carefully.
2. With clean hands, remove items from disinfectant solution, rinse under warm tap water, and air dry on a clean paper towel. Store in a zip-lock bag.

**NOTE**—Do not towel dry nebulizer parts; this could cause contamination.

### WARNING

To prevent possible risk of infection from contaminated cleaning solutions, always prepare fresh solution for each cleaning cycle and discard solution after each use.

3. Keep the outer surface of the tubing dust-free by wiping regularly. Nebulizer tubing does not have to be washed because only filtered air passes through it.

**NOTE**—The disposable nebulizer will last at least 15 days and possibly longer, depending on usage. Proper cleaning will also help extend the life of the nebulizer. Because it is disposable, we recommend that an extra nebulizer be kept on hand at all times. The reusable nebulizer is dishwasher safe and may be cleaned and reused up to 6 months.

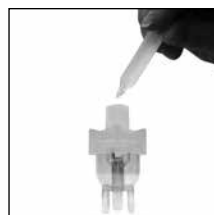
Figure 3

Figure 4

Figure 5

Figure 6

Figure 7

## COMPRESSOR CLEANING

---

1. With power switch in the "Off" position, unplug power cord from wall outlet.
2. Wipe outside of the compressor cabinet with a clean damp cloth every few days to keep dust-free.

### DANGER

Do not submerge in water; doing so will result in compressor damage.

## FILTER CHANGE

---

1. Filter should be changed every 6 months or sooner if filter discolors.
2. Remove filter (Fig. A) by grasping it firmly and pulling out of the unit. Discard old filter.
3. Replace with a new DeVilbiss filter (part #3655D-601). Additional filters should be purchased from your DeVilbiss provider.
4. Push new filter into position.

**CAUTION**– Do not reuse or wash filter or substitute any other material such as cotton for a DeVilbiss air-inlet filter. It will result in compressor damage.

## MAINTENANCE

---

All maintenance must be performed by a qualified DeVilbiss provider or authorized service center.

### DANGER

Electric shock hazard. Do not remove compressor cabinet. All disassembly and maintenance must be done by a qualified DeVilbiss provider.

## SPECIFICATIONS

---

Size	7.5" W x 4.0" H x 7.2" D (19.1 cm x 10.2 cm x 18.3 cm)
Weight	
3655D	3.8 lbs. (1.72 kg)
3655I	4.8 lbs. (2.18 kg)
Electrical Requirements	
3655D	2.5 amps, 115 VAC (±10%) 60 Hz
3655I	1.4 amps, 230 VAC (±10%) 50 Hz
Maximum Sound Pressure Level	56 dBA
Power Consumption	
3655D	140 watts maximum
3655I	160 watts maximum
Maximum Compressor Pressure	35 psig (241 kPa)*
Compressor Free Air Flow	13 lpm*
Operating Temperature Range	+40° to +104°F (+5 to +40°C)
Operating Humidity	up to 95% non-condensing
Storage/Transport Temperature Range	-40° to +158°F (-40° to +70°C)
Storage/Transport Humidity	up to 95% non-condensing
Supplied Nebulizer Capacity	10 ml (cc)
Supplied Nebulizer Nebulization Rate	0.15 ml/min or greater*
Supplied Nebulizer MMAD	5 microns or less**
Supplied Nebulizer Operating Pressure	12 psig (min)*

\*Conditions may vary based on altitude above sea level, barometric pressure, and temperature.

\*\*Value determined with minimum flow rate through the nebulizer (20°C, sea level).

There are no hazardous materials used in the construction of this unit, therefore, there are no known risks associated with disposal of the equipment.

This product (model 3655D) complies with the following electromagnetic compatibility standard: IEC 60601-1-2:2nd Edition

**3655D Equipment Classifications:**

Protection against electric shock: ..... Class II  
Degree of protection against electric shock: ..... Type BF  
Equipment not suitable for use in the presence of flammable anaesthetic mixture with air or with oxygen or nitrous oxide  
Mode of operation: ..... Continuous operation

**FIVE-YEAR LIMITED WARRANTY ON COMPRESSOR PARTS**

All DeVilbiss Compressor models are warranted to be free from defective workmanship and material for a period of five years from the date of purchase. Any defective part(s) will be repaired or replaced at DeVilbiss' option if the unit has not been tampered with or used improperly during that period. Make certain that any malfunction is not due to inadequate cleaning or failure to follow the instructions. If repair is necessary, contact your authorized DeVilbiss provider or DeVilbiss Service Department at (800-338-1988, 814-443-4881) for instructions.

**NOTE**— *This warranty does not cover providing a loaner compressor, compensating for costs incurred in compressor rental while said unit is under repair, or costs for labor incurred in repairing or replacing defective part(s).*

THERE IS NO OTHER EXPRESS WARRANTY. IMPLIED WARRANTIES, INCLUDING THOSE OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED TO THE DURATION OF THE EXPRESS LIMITED WARRANTY AND TO THE EXTENT PERMITTED BY LAW ANY AND ALL IMPLIED WARRANTIES ARE EXCLUDED. THIS IS THE EXCLUSIVE REMEDY AND LIABILITY FOR CONSEQUENTIAL AND INCIDENTAL DAMAGES UNDER ANY AND ALL WARRANTIES ARE EXCLUDED TO THE EXTENT EXCLUSION IS PERMITTED BY LAW. SOME STATES DO NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, OR THE LIMITATION OR EXCLUSION OF CONSEQUENTIAL OR INCIDENTAL DAMAGES, SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY TO YOU.

This warranty gives you specific legal rights, and you may also have other rights which vary from state to state.

**Manufacturer's Note**

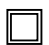
Thank you for choosing a DeVilbiss compressor/nebulizer. We want you to be a satisfied customer. If you have any questions or comments, please send them to our address on the back of this manual.

<p><b>For Service Call Your Authorized DeVilbiss Provider:</b></p> <p>_____</p> <p>_____</p> <p>_____</p> <p>Phone _____</p> <p>Purchase Date _____</p> <p>Serial # _____</p>
---

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DeVilbiss® and Pulmo-Aide® are registered trademarks of DeVilbiss Healthcare.
© 2021 DeVilbiss Healthcare LLC. 10.21 All Rights Reserved.
**Drive DeVilbiss Healthcare Ltd.**
Limited**
**DeVilbiss Healthcare GmbH**
DeVilbiss Healthcare LLC • 100 DeVilbiss Drive • Somerset, PA 15501 • USA