



Play. Engage.
Participate.



Inspired by Drive

Pediatric Product Line

www.inspiredbydrive.com





Table of Contents

Bathing

Marlin Bath Lift	1
Otter® Bathing Chair	2
Versa™ Transfer System	3
Ultima™ & Ultima Access Bath Chairs	4
Ultima™ & Ultima Access Transfer System	5-6

Car Seat

Spirit Plus™ Car Seat	7
Spirit™ Spica Car Seat	8

Wheeled Mobility

Miko® Tilt-In-Space Mobility Base	9
Trotter Mobility Chair	10
Wallaby Folding Pediatric Wheelchair	11

Walking Aids

Nimbo™ Posterior Walker	12
Universal Forearm Crutches	13
Moxie™ GT Gait Trainer	14

Seating & Positioning

P Pod®	15
Seat2Go™	16
First Class School Chair	17

Toileting

Contour™ Toilet Supports	18
Pinniped™ Pediatric Commode	19

Beds

The Bradshaw Junior™ Bed	20
--------------------------	----



Marlin Bath Lift

With postural support, the Marlin makes bath time enjoyable and safe for users with poor balance and trunk control who are unable to sit in a standard bathtub. The functional and easy-to-use bath lift features a reclining back and rechargeable battery-operated lift mechanism.

Features and Benefits

Safe Operation

- Waterproof, simple-to-use hand control safely and gradually lowers user into bathtub
- Safety mechanism prevents lift from lowering if battery charge is not sufficient to raise it back up

Adjusts to User Needs

- Seat height adjusts from 2.5–19" and back reclines from 10° off vertical to 50° at its lowest position
- Multiple positioning accessories available for postural support

Lightweight, Versatile & Travel Friendly

- Supports up to 300 lbs. and has a base unit weight of 20.5lbs
- Disassembles for easy storage and transport
- Tool-free assembly and installation

Accessories



Head and Trunk
Support Pads



Swivel and
Transfer Aid



Rapid Dry
Positioning Belt



Contoured
Chest Harness



The Otter® Bath Chair

The Otter® Bath Chair provides a safe, comfortable and adjustable bathing option to meet individual needs.

Features and Benefits

Versatile & Customizable Design

- Plastic fold-flat frame comes with standard or soft removable and washable fabric
- Multiple color and fabric options available to fit any bathroom décor
- Four frame sizes to choose from to meet user needs

Adjusts to User & Caregiver Needs

- Uni-bar allows for five one-handed angle adjustments while user is in or out of chair
- Trunk, leg positioning straps and lateral cushions provide head and trunk support
- Adjustable, slip-resistant legs or shower stand add extra height to accommodate bathing needs

Accessories



Lateral Supports




Shower Stand



Tub Stand




Standard Fabric

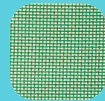

Vinyl
Seahorse
Green

Soft Fabric

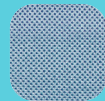

Beach
Bubble
Blue



Lionfish
Lavender



Versa™ Transfer System

The lightweight and multifunctional Versa™ Transfer System is most suitable for users who require minimal positioning support. Constructed of high-grade aluminum, stainless steel and padded polyurethane, the Versa can be utilized as a bath transfer system, roll-in shower chair or bedside commode.

Features and Benefits

Versatile Design

- Compact 18" x 18" tub base transfer frame fits inside most bathtubs
- Multiple bridge options to fit your bathroom: Standard removable 13" bridge, Optional fixed bridge in 9", Optional adjustable bridges from 13½–17¾" or 17½–21¾" lengths

Multifunctional & Built-In Safety

- Allows incontinent users to bathe and shower comfortably with an opening for easy toileting or bathing access
- Locking system limits movement and slipping

User Accommodation

- Use for adolescents, teens or adults up to 300 lbs.

Accessories



Fixed-Length
Bridge



Adjustable
Length Bridge



Rapid Dry
Positioning Belt



Lap Belt



Ultima™ & Ultima Access Bath Chair

Designed with comfort, functionality and durability in mind, the Ultima™ and Ultima Access Bath Chairs provide great solutions for those with assistive bathing needs. The Ultima and Ultima Access feature a safe bathing chair with independently adjustable angles, available in multiple sizes, two colors and a durable metal frame.

Features and Benefits

Multiple Sizes to Choose From

- 3 sizes available in the Ultima, 2 sizes in the Ultima Access
- Ultima Access has ease of access opening for cleaning and/or toileting
- 3 weight caps: 130lb/160lb/180lb

Easy to Keep Clean

- Bacterial and fungal-resistant ComfortMesh™

Improved Patient Comfort and Safety

- 180° back recline
- Multiple support pads and positioning belts available

Adjustable Seating

- Independently adjustable seat, back and leg angles
- Frame folds for easy storage

Accessories



Detachable
Extension Base



Detachable
Rolling Base



Positioning
Belt



Rapid Dry
Positioning Belt



Beach
Bubble
Blue



Lionfish
Lavender



Ask about our Contour
Deluxe Bath Chair



Ultima™ & Ultima Access Bath Transfer System

The Ultima™ / Ultima Access Bath Transfer Systems are designed to significantly reduce the effort required by the parent or caregiver to transfer a child or adult into a tub for bathing.

Features and Benefits

Improved Patient Comfort & Safety

- Multiple adjustment angles for comfort
- Comes standard with ComfortMesh™ seat, tub base, rolling base, two positioning belts and anti-tippers.
- Available in two sizes: S/M or L/XL
- Ultima Access has ease of access opening for cleaning and/or toileting
- Many optional comfort and positioning accessories available.

Available Colors



Beach
Bubble
Blue



Lionfish
Lavender

Accessories



Footrests



Attached
Rolling Base



Foldable
Transfer Base



Compact
Transfer Base



Head Support

Functional & Durable

- Tub bases are constructed with high-grade aluminum and stainless steel which provides exceptional strength and durability.
- The Slider Safe™ Auto-Stop Safety System ensures that the chair will not come off the base.
- Bacterial and fungal resistant ComfortMesh™ seat for easier cleaning


Two Tub Base Options: Folding or Compact

- Folding Base places two legs in the tub and two outside the tub, and folds for easy storage.
- Compact Base places all four legs inside the tub for use in tighter spaces, and is non-foldable.

Easy to Use

- Three components: a mesh seating system for comfort, a rolling base used for the initial portion of the transfer, and a tub base that when connected to the rolling base engages the shuttle function for transfer.
- Connecting the rolling base to the tub base is easy. Simply line up the interlocking tubes and secure a latch. No tools required.

Multifunctional bath transfer system for use in home or facility without any bathroom modifications. All systems include anti-tippers.

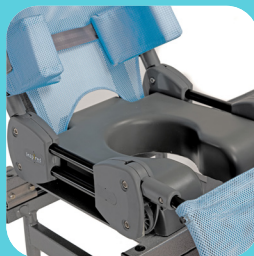

Ultima Access model comes with polyurethane seating area which creates an opening in the chair allowing the system to function as a commode.



Lateral Support
Pads



Rapid-Dry Pelvic
Positioning Belt



Pelvic
Positioning
Belt



Detachable
Armrest



Ask about our
Spirit Standard

Spirit Plus™ Car Seat

Provides positioning support for users with mild to complex positioning needs. The Spirit Plus™ is the only car seat on the market today offering swing-away padded trunk supports, hip supports and abductor.



Features and Benefits

User Accommodation

- Two sizes with seat depth extender (12"–16") to accommodate growth
- Height & weight capacities support users up to 66" and 130 lbs.

Optimal Support

- Recline bar adds 25° of extra tilt for added comfort
- Trunk & hip supports and abductor help maintain proper positioning postural support

User Safety & Security

- Memory foam back/seat padding reduces sores and fatigue
- Open shell design with low-profile sides and swing-away supports for easy transfers

Available Colors



Convertible
Pink



Go Kart
Green



Rally
Blue



Roadster
Red



Speedway
Gray



Incontinence
Covers



3-Inch Extensor
Thrust Wedge



Buckle Guard



Retainer Clip
Guard



Spirit™ Spica Car Seat

The Inspired by Drive Spirit™ Spica Car Seat features a low and wide profile with positioning wedges to safely transport users with spica casts. The Spirit Spica seat with wipeable cover accommodates a wide range of user heights and weights to meet individual needs.



Features and Benefits

Built-In Safety

- 5-point adjustable seat harness with included tether strap
- Easily installed in most vehicles
- Two seat depth options with a wide weight and height range

Adapts to User

- Seat depth options of 12" and 16"; accommodates users up to 66" tall and a 25–130 lb weight capacity
- Adjustability and seat configurations can be obtained for different sized clients and Spica casts

Comfort & Convenience

- High-quality memory foam on the back and seating surface provides a comfortable environment

Included with the Spirit™ Spica Car Seat



Positioning
Wedges (2)



Memory foam
Support



Wipeable, low-maintenance
upholstery



Miko® Tilt-In-Space Mobility Base

Available in 3 fun colors, the Miko offers a sporty frame, stylish flexible seating, and superb positioning support for children with limited mobility.

Features and Benefits

Lightweight & Durable

- Frame folds for transport
- 2 sizes available, seating system accommodates growth to help maximize results

Improved Patient Comfort & Safety

- Easy to adjust, 40 degrees of posterior tilt and 5 degrees of anterior tilt
- 180 degrees of back recline

Crash-Tested for Transportation

- Ensures safe transport

Available Colors



Fire Truck
Red



Jet Fighter
Blue



Punch Buggy
Pink

Accessories



All-Terrain
Option



Oxygen
Tank
Holder



Angle-
Adjustable
Footbox



Upper
Extremity
Support



Poziform
Belts &
Harnesses

Trotter® Mobility Chair

The Trotter® Mobility Chair offers parents and caregivers a lightweight, foldable alternative to heavier wheelchair-style mobility bases.

Features and Benefits

Improved User Comfort & Safety

- Poziform® 5-Point Harness with padded breathable covers provides additional support
- Crash tested for transport to WC 19 Standards
- Adjustable seat depth, back-to-seat angle and tension positioning straps for user customization and growth

Easy Maneuverability & Transport

- Curb assist pedal, 8" front casters and 10" rear wheels help navigate curbs and all types of terrain
- Lightweight and foldable frame with closure strap for easy transport
- Canopy and under-seat carrier included

Accessories



Headrest
Extension



Contoured
Wings



Lateral
Support
& Scool Strap



Upper
Extremity
Support



Full
Torso
Vest



Rear
Anti-Tippers



Available Colors



Fire Truck
Red



Jet Fighter
Blue



Punch Buggy
Pink



Wallaby Folding Pediatric Wheelchair

Perfect for active users, the Wallaby™ Foldable Pediatric Wheelchair provides mobility, comfort and safety.

Features and Benefits

User Customization

- 12" Width Red and 14" Width Blue

Safe & Easy Transport

- WC-19 Compliant when used with Headrest Extension
- Folds for storage

Available Colors



Fire Truck Red



Jet Fighter Blue



Also available the
Tilt-In-Space Kanga
Wheelchair

Accessories



Headrest Extension



Chest Harness



Ankle Positioners



Elevating Legrests





Nimbo™ Posterior Walker

The Nimbo™ Posterior Walker provides assistance and postural support for safe ambulation and is lightweight, height-adjustable and foldable.



Features and Benefits

Built-In Growth

- Lightweight and durable aluminum frame folds easily for transport and storage
- Five frame sizes with 4" or more of adjustment for a range of growth in each size

Safe & Smooth Mobility

- 5" rubber front swivel wheels can be locked in a forward position
- Large wheels provide smooth navigation indoors and outdoors
- Metal, anti-rollback assembly provides improved durability

Flexible Customization

- Versatile design allows for anterior configuration if needed

Accessories



Soft Seat
Harness



Pelvic
Stabilizer



Forearm
Platforms



Fold-Down
Seat



9" Front
Wheels



Knight Blue



Castle Red



Wizard Purple



Emperor Black
(Extra-Large
Only)



Universal Forearm Crutch

The lightweight Inspired by Drive Universal Forearm Crutches provide a comfortable alternative to standard crutches and canes for users who need moderate support while standing.

Features and Benefits

Durable & Adjustable

- Aluminum construction with push-button height adjustment in 1-inch increments
- Full-style composite forearm cuffs adjust individually to accommodate each user
- 1-year limited warranty

Available Colors



Knight
Blue



Wizard
Purple



Moxie™ GT Gait Trainer

The Moxie™ GT offers a tremendous amount of versatility to meet the needs of most users. Height-adjustable and available in three sizes, the Moxie GT can be used in either a posterior or anterior position.



Features and Benefits

Customizable Frame Options

- Lightweight and durable aluminum and steel frame construction
- Medium and large frames fold easily for storage and transportation
- Three frame sizes with 1" height adjustments allow for range of growth in each size

Safe Ambulation

- Large 5.5" rubber multi-function front and rear casters smoothly navigate indoors and outdoors
- Wheel locks, swivel locks, unidirectional controls and drag control for added safety
- Accommodates anterior or posterior configurations

Accessories



Hand Loops



Forearm Platforms



Handgrips



Guide Bar



Trunk Support



Warrior Blue
(Small Only)



Magic Green
(Medium Only)



Sword Gray
(Large Only)



P POD®

The P Pod® provides clinical benefits to the user while encouraging active participation at home or school. The specially molded foam liner provides outstanding positional support and comfort. The P Pod consists of the positioning pod, molded foam liner for optimal postural support, bolster support for upright positioning and bean bag base.



Features and Benefits

Support & Positioning

- Provides midline alignment of trunk, pelvis and lower extremities
- Molded foam liner offers clinically appropriate postural support
- Pelvic belt keeps user in desired position; headrest provides added head positioning

Lightweight & Customizable

- Adjustable recline angle based on pod alignment on bean bag base
- Lightweight design features removable and washable fabric cover or Dartex™ wipeable fabric
- 180° of back recline

Accessories



Activity
Frame



Chest
Harness



Mobility
Base



Foot & Lateral
Bolster



Lap/Freestanding
Tray

Available Colors



Mermaid
Pink



Pool Blue
(Small & Medium
Only)



Whale Blue
(Large & Extra-
Large Only)



Jellyfish Purple
(Large & Extra-
Large Only)



Starfish Red
Dartex



Seat2Go™

The lightweight and portable Seat2Go™ positioning seat provides postural support and can be secured safely and easily to any chair or stroller.



Features and Benefits

Safe & Comfortable

- Contoured seat prevents user from sliding forward and improves seating posture
- Double-strap trunk support and positioning belt included
- 95° seat-to-back angle
- Two sets of straps attach to any chair
- Seat straps store under seat behind clear flap and back straps store in back pocket

Multi-Use Design

- Can be used as an insert for a single or double stroller or as a dining/activity chair
- Lightweight and versatile

Accessories



Height-Adjustable Headrest



Depth-Adjustable Abductor



First Class School Chair

The First Class School Chair is designed to grow with your child and is available with multiple positioning accessories to help users stay safe while seated in the classroom, home or during therapy. A specially designed chair can provide safety and good posture unavailable with standard school and office chairs.

Features and Benefits

Long-Term Durability

- Metal construction with high-grade vinyl and foam
- 1-year limited warranty

Safe & Supportive

- Seat belt provides pelvic support and maintains proper hip positioning
- Positioning equipment provides an extra level of safety while sitting
- Standard padded hip belt included

Adjustable

- Two adjustable frame sizes to choose from
- Height / depth-adjustable seat and adjustable-height flip-up armrests

Accessories



Headrest
Extension



Chest
Harness



Ankle
Positioners



Elevating
Legrests



Anti-Tippers



Contour™ Toilet Supports

The Contour™ Series Toilet Supports Series is designed to assist with positioning and comfort through the use of support belts and padding.



Features and Benefits

Built-In Growth & Comfort

- Padded back reduces pressure
- Orthopedically contoured back supports can be added for increased comfort and support
- U-Shape Mounting Base allows you to upgrade or replace the support back to the next size for a fraction of the cost
- Lo-Back and Hi-Back options available

Durable Construction

- Stainless steel heavy duty mounting base
- Fast and easy to clean with a wipeable seat
- Latex and silicone-free

Accessories



Armrest & Foot-rest



Padded Reducer Ring



Plastic Reducer Ring



Soft-Flex Splash Guard



Pinniped™ Pediatric Commode

The Pinniped™ Pediatric Commode includes an all-in-one welded steel frame with a plastic seat and contoured back for added comfort and durability.

Features and Benefits

Adapts to User Needs

- Multiple height adjustments allow user to place feet on floor for comfort, safety and security
- Adjustable hook and loop chest strap for added support
- Rubber non-skid tips for floor protection

Durable Comfort

- Molded plastic armrests and contoured back
- Snap-on plastic seat

Hygienic Solutions

- 7.5 qt. removable commode bucket with metal handle and cover
- Optional commode pail liners help prevent splashing and keep bucket clean

Accessories



Biodegradable
Commode Liners



Poziflex® Chest
Straps



Bradshaw Junior™ Fully Electric, Height-Adjustable Pediatric Bed

Safer Sleep Positioning for Toddlers to Teens

Designed for toddlers to teens and everyone in between, the Bradshaw™ Junior brings you clinical positioning without the clinical look. The fully electric design offers effective sleep positioning for children and young adults up to 363 lbs., including features like foot elevation and Reverse Trendelenburg positioning, and an adjustable height range allows comfortable positioning for caregivers.



Setup is simple: each section weighs under 44 pounds and rolls easily on a dedicated transport stand, so one person can assemble it with ease. Optional double-lock side rails add peace of mind, and a handset lockbox keeps controls secure.

Features and Benefits

Comfortable Positioning

- Offers foot elevation and Reverse Trendelenburg tilt to support circulation, comfort, and effective care
- Variable height range helps caregivers assist safely and comfortably.

Home-Friendly Design

- Finished in a clean oak wood tone that blends seamlessly into bedrooms

Easy to Setup & Move

- Includes a rolling transport stand for easy moving
- Each section weighs under 44 pounds for easy assembly by one person

Safety & Security

- Optional full-length side rails include a two-stage lock for added protection
- Lockbox keeps controls safe from tampering



877-224-0946 | www.inspiredbydrive.com

©2025 Medical Depot, Inc. dba Drive DeVilbiss Healthcare

CA_rev_09_12_2025